



楚天微球
TRUKING MICRO-SPHERE

**Truking Micro-sphere Biotechnology Co.
Product manual**

Instruction Manual No.: 610
Edition number: 01
Effective date: 2025.01.01

TK-Col 6FF

Gel Filtration Prepacked Columns


Product Manual



楚天微球生物技术(长沙)有限公司
TRUKING MICRO-SPHERE BIOTECHNOLOGY (CHANGSHA) CO., LTD

Content

1. Product Introduction	3
2. Technical parameters	3
3. Methods of use	6
4. Cleaning and regeneration	8
5. Sterilization and storage	8
6. Destruction and recycling	9
7. Ordering Information	9

 楚天微球 <small>TRUKING MICRO-SPHERE</small>	Truking Micro-sphere Biotechnology Co. Product manual	Instruction Manual No.: 610 Edition number: 01 Effective date: 2025.01.01
---	--	---

1. Product Introduction

The TK-Col 6FF series of pre-packed gel filtration columns feature TA-6FF agarose gel filtration media packed into TK-EC chromatography columns (TK-EC 16/70, TK-EC 16/100, TK-EC 26/70, TK-EC 26/100) into ready-to-use gel filtration columns. This eliminates the hassle of self-packing columns and the risk of suboptimal column efficiency. The perfect integration of media and column makes chromatography efficient and straightforward. The TK-Col 6FF series gel filtration pre-packed columns are suitable for industrial process development and large-scale sample preparation, ideal for the separation and purification of macromolecules such as viruses, polysaccharides, and nucleic acids.

2. Technical parameters

Table 1 Technical parameters

Product name		TA- 6FF
Agarose concentration		6%
Separation range	Linear molecules	10kd~1000kd
	Spherical molecules	10kd~4000kd
	Nucleic acid	15bp~190bp
Particle size distribution range		45-165μm
Average particle size		90μm
Recommended Flow Rate		110cm/h (TK-EC50/100, h=80cm)
Maximum Pressure		0.3MPa (3bar)
pH Stability		2-12 (working),2-14 (CIP, short-term)
Chemical stability		Stable in common aqueous solutions: 2M NaOH, 8M urea, 6M guanidine hydrochloride, 1M acetic acid, 70% ethanol, 30% isopropanol
Sterilization		Resistant to 121° C, 20min sterilization
Storage conditions		20% ethanol, room temperature (2~8° C is better)
Appearance		White paste, can be layered when placed
Precautions		Freezing destroys the internal structure; 2% Benzyl Alcohol for foreign shipments only or customer specified

Table 2: Technical parameters for each pre-assembled column (see end page for article number)

Product name	TK-Col 16/60 6FF	TK-Col 16/90 6FF	TK-Col 26/60 6FF	TK-Col 26/90 6FF
Item No.	Y60101	Y60102	Y60103	Y60104
Resin	TA-6FF	TA-6FF	TA-6FF	TA-6FF
Separation range (Da)	10-4000	10-4000	10-4000	10-4000
Average particle size	90μm	90μm	90μm	90μm
Column bed volume (ml)	120	200	320	530
Maximum Pressure	0.3	0.3	0.3	0.3
Recommended flow rate (ml/min)	4-13	4-13	10-35	10-35
Storage	2-30℃ , 20% ethanol or 2% benzyl alcohol (for international transportation)			

3. Methods of use

- ◆ The chromatography column is made of glass and should be handled gently to prevent it from breaking or affecting the column efficiency.
- ◆ To avoid clogging the column, all samples and buffers need to be filtered with 0.45um membrane.
- ◆ In order to get a good separation effect, avoid too much temperature change between the buffer and the chromatography column.
- ◆ Place the chromatography column in a place without direct sunlight.
- ◆ The chromatography column can be used in a chromatography cooler, but the flow rate needs to be reduced appropriately.

3.1 Connecting the column to the chromatography system

- Open the package and take out the chromatography column.
- Check whether the chromatography column is intact, and whether the chromatography column has been dried out by air intake during transportation, if any of the above situations occurs, please contact the sales representative of Truking Micro-sphere in time.
- Fix the chromatography column next to the chromatography system, pay attention to the flow direction of the chromatography column.
- Start the chromatography system, make sure the air bubbles in the chromatography system are drained, and set the alarm pressure of the chromatography system to 0.3MPa, then adjust and keep the flow rate running at 0.2ml/min.
- After the chromatography system is purged of air bubbles, open the upper and lower plugs of the chromatography column and connect the chromatography column under low flow rate operation.

3.2 Pretreatment of chromatographic columns

- Rinse, the chromatography column is stored in 20% ethanol or 2% benzyl alcohol (for international shipments) during transportation, first rinse off the storage solution with 2 column volumes of distilled water.
- Sterilization, for sample safety, it is recommended to rinse 2 column volumes with 0.5M NaOH and then 2 column volumes with distilled water for the first use.

3.3 Balancing

- The recommended flow rate is generally used to rinse the column with an equilibration buffer. The choice of equilibration buffer depends on the stability of the sample, and the type and pH of the buffer do not have a significant effect on the effectiveness of gel filtration. However, agarose contains a small amount of sulfate and carboxylate groups and in order to minimize non-specific adsorption of basic protein samples, it is recommended that at least 0.15 M NaCl be added to the equilibration buffer. The pH and conductivity of the buffer to be exported should be the same as that of the buffer before entering the column to indicate that the column is well equilibrated, usually 2-5 column volumes are required.

- Use the recommended flow rate to flush the chromatographic column with equilibrium buffer. The pH and conductivity of the buffer to be exported are consistent with that of the buffer before entering the chromatographic column that means the chromatographic column is well equilibrated, generally 2~5 column volumes are needed.
- Depending on the type of column, the recommended flow rates are generally chosen within the range of flow rates, see Table 2.

3.4 Sampling

- In order to prevent the sample from clogging the column, the sample needs to be filtered with 0.45μm microporous filter membrane before loading, and can be appropriately diluted if the viscosity is too high, and the protein concentration should not exceed 70mg/ml.
- Load the sample onto the chromatography column through the sample loading ring of the chromatography system and other devices, the volume of the sample will vary according to the difference between the target material and the impurity. If the difference between the target and the impurity is relatively small, the sample volume is 1-5% of the column volume, and if the difference between the target and the impurity to be removed is relatively large (e.g., removal of inorganic salts or buffer exchange in the sample), the sample volume is 10-25% of the column volume.

3.5 Separation

- Continue to flush the chromatography column with equilibration buffer, collecting the different components of the effluent based on UV absorption, until no more biomolecules are effluent, typically requiring 1 to 2 column volumes.

3.6 Rebalancing

- Rinse the column with equilibrium buffer and wait until the pH and conductivity of the outlet buffer are the same as the equilibrium buffer prior to entering the column before proceeding to the second sample, and so on.

3.7 Regeneration


- Flush the column with a solution containing 1M NaCl for 1 to 2 column volumes to rinse out nonspecifically adsorbed impurities.

3.8 Column Effectiveness Evaluation

Column efficiency can be determined by using acetone as indicator or NaCl as indicator, and the indicator solution and mobile phase are prepared according to the following table.

Table 3: Column efficiency determination methods

Methods	Acetone Method for Column Efficacy	Column Efficacy by NaCl Method
Sample	1.0% (v/v) acetone in water	0.8M NaCl (dissolved in water)
Sample volume	1.0% column volume	1.0% column volume
Mobile phase	Water	0.4M NaCl aqueous solution
Flow rate	30 cm/h	30 cm/h

 楚天微球 TRUKING MICRO-SPHERE	TruKing Micro-sphere Biotechnology Co. Product manual	Instruction Manual No.: 610 Edition number: 01 Effective date: 2025.01.01
Detection Data	UV 280 nm	Conductivity

3.9 Calculating Column Effect

Theoretical plate height (HETP), theoretical number of plates (N) and asymmetry factor (As) were calculated from the UV or conductivity curves with the following equations:

$$HETP = L/N$$

$$N = 5.54(V_R/W_h)^2$$

Where: V_R = retained volume

W_h = half peak width

L = column height


N = theoretical plate number

The units of V_R and W_h should be the same;

$$A_s = b/a$$

Where: a = first half peak width at 10% peak height

b = second half peak width at 10% peak height

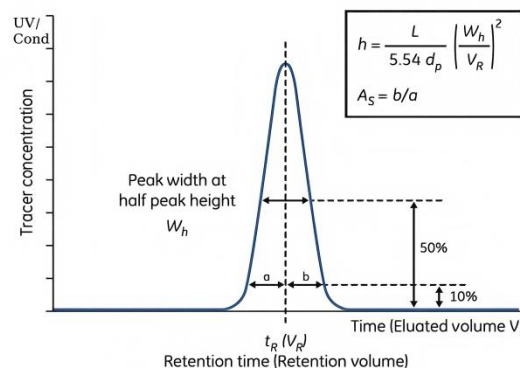

4. Cleaning and regeneration

TK-Col series gel filtration preloaded columns in use for a period of time there may be a decline in column efficiency, the separation effect deteriorates, the need for cleaning and regeneration, generally 5 ~ 10 cycles need to do a thorough regeneration (regeneration frequency depends on the contamination status of the columns), regeneration needs to be based on the nature of the contaminants using the appropriate regeneration reagents.

- First flush 1 column volume with buffer containing 1M NaCl
- To remove denatured proteins: 2 column volumes were backflushed with 1M NaOH at a flow rate of 20cm/h.

Note : Denatured proteins can also be removed by protease using 1mg/ml of gastric enzyme dissolved in 0.1M acetic acid solution containing 0.5M NaCl.

- Removal of lipids or lipoproteins: 70% ethanol or 30% isopropanol at a flow rate of 40 cm/h for 4 column volumes (to prevent air bubbles can be used to gradually increase the proportion of organic solvents in a gradient); or 1% non-ionic decontaminants
- Inorganic contaminants: rinse 2 column volumes with 0.5M acetic acid
- Rinse 4 column volumes with distilled water

5. Sterilization and storage

Since 20% ethanol or 2% benzyl alcohol preservation solution can inhibit the growth of bacteria, but does not have the effect of sterilization and pyrogen removal, it is recommended that TA-6FF media can be washed with 0.5M NaOH at the recommended flow rate on the loaded columns prior to and during use to reduce the risk of microbial contamination.

TA-6FF media is sold with 20% ethanol or 2% benzyl alcohol as preservation solution.

TA-6FF media are stored in 20% ethanol in a closed container at 2-30°C. To prevent ethanol evaporation and

microbial growth, it is recommended that the preservative solution be replaced with a fresh one every 3 months.

6. Destruction and recycling

TA-6FF is difficult to degrade in nature and incineration is recommended to protect the environment.

7. Ordering Information

Table 5 Article number and packaging

Product	Item No.	Norm
TK-Col 16/60 6FF	Y60101	1 piece
TK-Col 16/90 6FF	Y60102	1 piece
TK-Col 26/60 6FF	Y60103	1 piece
TK-Col 26/90 6FF	Y60104	1 piece